



**MEDICATECH USA®**

MedicatechUSA presents  
**EXPRIMER 17X17 Tethered DR Panel**

- ✓ Increased Dose Efficiency and Image Quality
- ✓ Improved Patient Throughput and Fast Workflow
- ✓ High DQE / MTF Value
- ✓ TRUVIEW®ART, Advanced Reverse Filtering Technology
- ✓ User Friendly & Ergonomic Design



**Single  
LAN Cable**




**Fast  
Workflow**



**Handle &  
Bracket**




**Single  
LAN Cable**

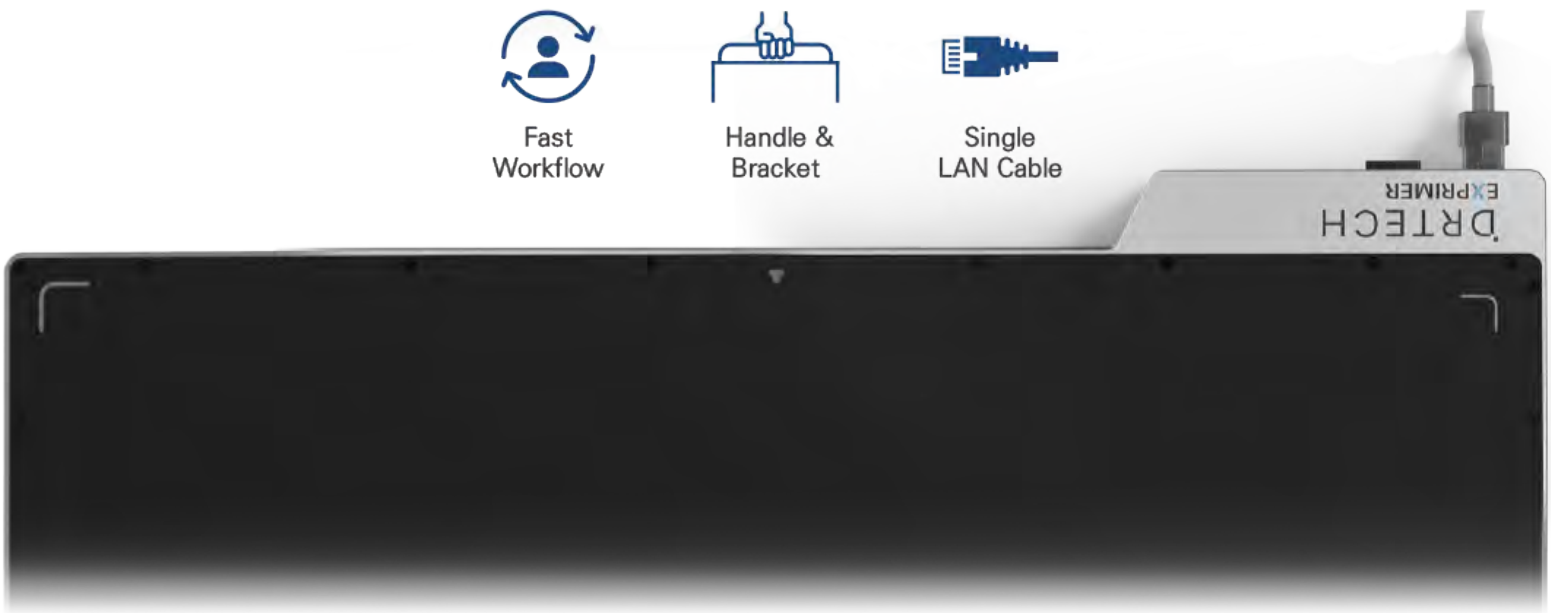



### Excellent Dose Reduction

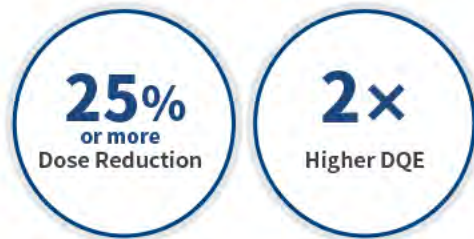


### Exceptional Diagnostic Clarity



### Focus Scintillator Technology

FOCUS (Fluorescent Optical Csi Upgraded Structure) has more focused optical signal compared to conventional Csi. As a result, ultra low noise is achieved by FOCUS technology.

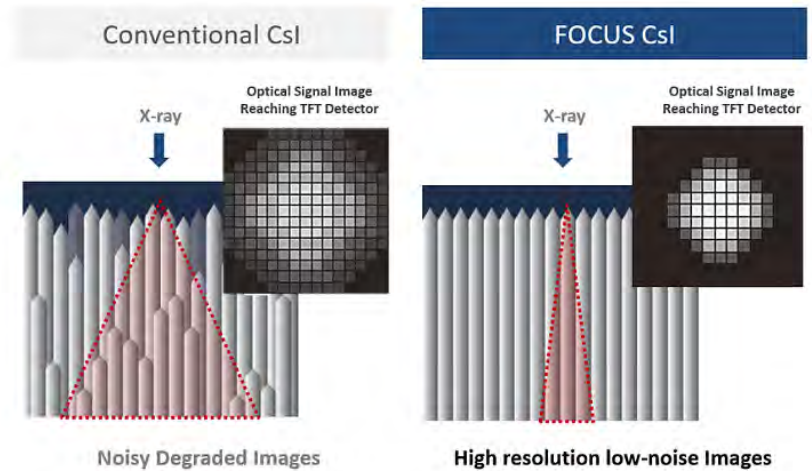

Excellent Image Capturing capacity for all radiographic procedures and packed in a well-designed and easy to use casing for many years of use. Ideal for any radiographic needs.



## Fast & Efficient Workflow


EXPRIMER Advanced series allows fast and efficient workflow brought by its fast exam time, speeding workflow and decreasing repeats.

First Image: Min. 4 sec.



## Rugged Design with Portability

EXPRIMER Advanced series' LAN cable with fixed protector increases the reliability and durability. Moreover, attachable handle allows carrying detector more conveniently.

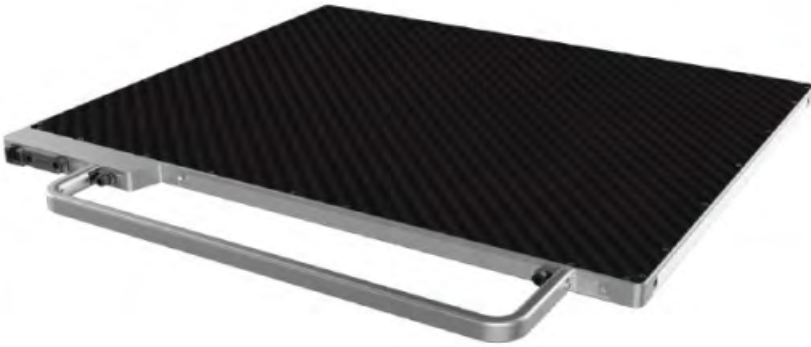

Model	<b>PL-DR-102</b>
Csl / Gadox	Csl / Gadox
Active area	430 x 430 (mm) / 17 x 17 (inch)
Resolution	3,072 X 3,072
Pixel pitch	140 $\mu$ m
Dimension	460 X 460 X 15 mm (Protruding Area: 24.2 x 125.5mm)
Weight	4.3 kg
X-ray Sync	Lossless AED/AWC, Sync Trigger
Data Interface	Gigabit Ethernet (LAN Cable)
ADC	16-bit

# VOYANCE®

By Mediatech USA®

An industry leading complete software suite that delivers a secure and cost-effective way to acquire, view and share DICOM Images. Voyance® is one of the easiest, most user friendly interfaces. Its modern design & attention to the user experience along with its and features & functions make it a great tool to every user.

**Instant installation and auto-configuration**

**Easy, user-friendly Modern user interface**

**Interface consistency across acquisition and viewing platforms**

**Optimum Workflow designed for efficiency**

**ClairVoyance® Image Processing Algorithm**

**Supports a wide range of detectors and generators**

**X-RAY Generator integration (CPI & DRGEM)**

**Auto Stitching and Auto Cropping**

**Side by side comparisons**

**Fully DICOM 3.0 Compatible: MWL, MPPS, and all DICOM ops.**

**Reject Analysis**





**MEDICATECH USA®**



**T. +1 800-817-5030**

**T. +1 949-679-2881**

**F. +1 949-679-2882**



**[www.medicattechusa.com](http://www.medicattechusa.com)**



**[info@medicattechusa.com](mailto:info@medicattechusa.com)**